

TRAPPER™

Exchange Device



Treat

efficiently in a 6–8F* guide
for maximum utility

Trap

confidently with a high-pressure
device that's recognizable at a glance

Repeat

multiple times leveraging a
tip-to-hub core wire for durability


TRAPPER Exchange Device is the only dedicated product engineered to allow treatment, then facilitate device exchanges without losing wire position to maximize the efficiency of your complex cases.



Master the Complex™

TRAPPER™

Exchange Device

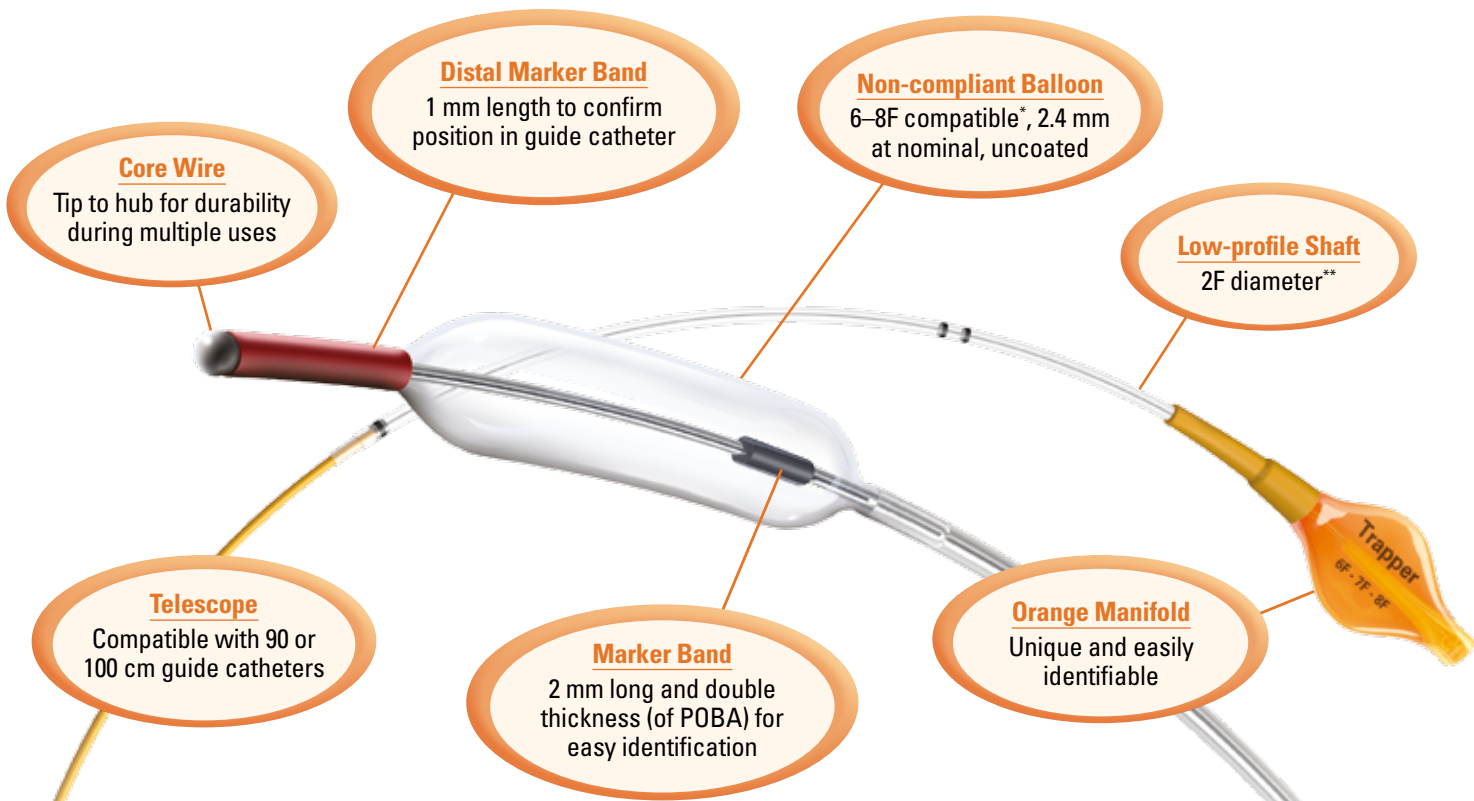

A fixed wire catheter that has a 10 mm balloon near the distal tip, designed to secure a 0.014" (0.36 mm) wire within a 6F through 8F guide catheter when inflated to 12 ATM.

There is a 1 mm marker band at the tip and a 2 mm marker band at the proximal balloon edge to facilitate **easy** placement and identification during **rapid** exchanges.

* 6F – 2.00 mm, 7F – 2.33 mm, 8F – 2.67 mm

** 2F – 0.67 mm

Balloon Length: 10 mm
Nominal Pressure: 12 ATM (1216 kPa)
Rated Burst Pressure (RBP): 20 ATM (2027 kPa)

Description	Order Number (GTIN)	Ref/Catalog Number
Trapper Exchange Device	05391526210154	H749 3933013 0

**Boston
Scientific**

Advancing science for life™

Caution: The law restricts these devices to sale by or on the order of a physician. Indications, contraindications, warnings, and instructions for use can be found in the product labelling supplied with each device. Information for the use only in countries with applicable health authority product registrations. Information contained herein is for use or distribution outside the US, France, and Japan only. All trademarks are the property of their respective owners.

IC-466604-AA APR 2017

www.bostonscientific.eu

© 2017 Boston Scientific Corporation or its affiliates. All rights reserved.
DINCAR2709EA